



SOMATEX[®]
Medical Technologies GmbH

Biopsy Handy

Semi-Automatic Biopsy Device

S A F E T Y T O T H E M A ^x

Biopsie-Handy Spezifikationen

Das SOMATEX® **Biopsie-Handy** arbeitet halb-automatisch und dient zur Gewinnung von histologisch verwertbarem Gewebematerial aus unterschiedlichen Weichteilen und Organen. Mit der **Biopsie-Schleuse** ist eine mehrfache Gewebeeinnahme mit demselben Instrument möglich.


Speziell für die Anwendung unter MR-Kontrolle gibt es das Biopsie-Handy auch als **MRT-System**.

Vorteile im Überblick

- ▶ Bequeme Einhandbedienung
- ▶ Geringes Gewicht
- ▶ Ultraschall- und CT-gestützt einsetzbar
- ▶ Als MRT-System unter MR-Kontrolle einsetzbar
- ▶ Atraumatische Nadelspitze
- ▶ cm-Markierung
- ▶ Optional mit koaxialer Schleuse erhältlich




Kanüle mit cm-Markierung
Cannula with cm-markers

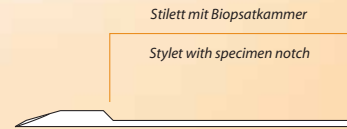

Stilett mit Biopsatkammer
Styilet with specimen notch

BIOPSIE-Schleuse
BIOPSY Sheath

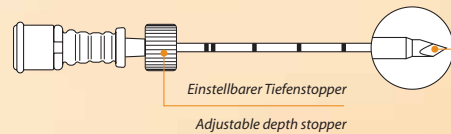

Trokar-Schliff
Trocar tip
Einstellbarer Tiefenstopper
Adjustable depth stopper

Biopsie-Handy / Biopsy Handy			
Artikelnummer Order no.	Gauge	Durchmesser Diameter	Länge Length
900 118	20	0,90 mm	100 mm
900 120	20	0,90 mm	150 mm
900 122	20	0,90 mm	200 mm
900 126	18	1,20 mm	100 mm
900 128	18	1,20 mm	150 mm
900 130	18	1,20 mm	200 mm
900 134	16	1,60 mm	100 mm
900 136	16	1,60 mm	150 mm
900 138	16	1,60 mm	200 mm
900 142	14	2,10 mm	100 mm
900 144	14	2,10 mm	150 mm
900 146	14	2,10 mm	200 mm

Biopsie-Schleuse / Biopsy Sheath			
Artikelnummer Order no.	Gauge	Durchmesser Diameter	Länge Length
170 251	18	1,15 mm	50 mm
170 253	18	1,15 mm	100 mm
170 255	18	1,15 mm	150 mm
170 257	17	1,50 mm	50 mm
170 259	17	1,50 mm	100 mm
170 261	17	1,50 mm	150 mm
170 263	15	1,85 mm	50 mm
170 265	15	1,85 mm	100 mm
170 267	15	1,85 mm	150 mm
170 269	13	2,25 mm	50 mm
170 271	13	2,25 mm	100 mm
170 273	13	2,25 mm	150 mm


Passende Schleuse
Fitting sheath

The SOMATEX® **Biopsy Handy** is a semi-automatic tool for obtaining histologically usable tissue material from a variety of soft tissues and organs. When complemented with the **Biopsy Sheath**, one Biopsy Handy can be reused for multiple tissue sampling.

Especially for application with MRI Control is the **MRI System** available.

Advantages at a glance

- ▶ Easy handling with one hand
- ▶ Light weight
- ▶ Ultrasound and CT compatible
- ▶ MRI System for application with MRI Control
- ▶ Atraumatic needle tip
- ▶ cm-markers
- ▶ Optionally available with coaxial sheath

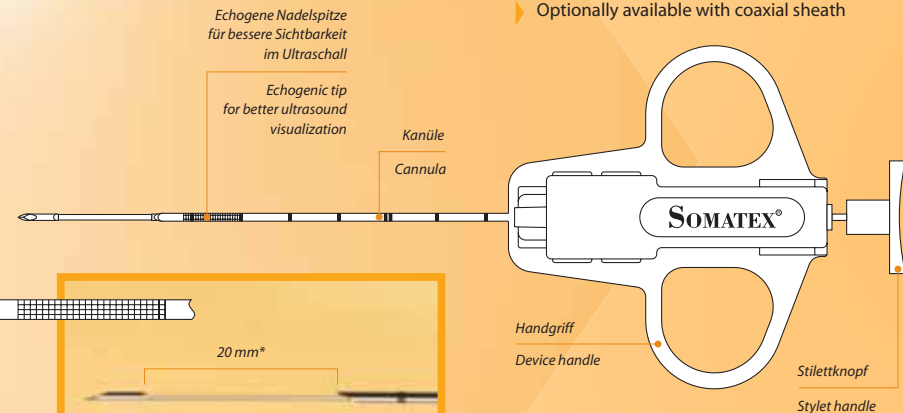

Echogene Nadelspitze
für bessere Sichtbarkeit
im Ultraschall

Echogenic tip
for better ultrasound
visualization

Kanüle
Cannula

Handgriff
Device handle

Stilettknopf
Styilet handle

*bei MRT-Biopsie-Handy abhängig vom Durchmesser
*with MRI Biopsy Handy dependent of diameter

MRT-Biopsie-Handy / MRI Biopsy Handy				
Artikelnummer Order no.	Gauge	Durchmesser Diameter	Länge Length	Biopsatkammer Specimen notch
603 108	18	1,20 mm	100 mm	10 mm
603 110	18	1,20 mm	150 mm	10 mm
603 112	18	1,20 mm	200 mm	10 mm
603 118	16	1,60 mm	100 mm	15 mm
603 120	16	1,60 mm	150 mm	15 mm
603 122	16	1,60 mm	200 mm	15 mm
603 128	14	2,10 mm	100 mm	20 mm
603 130	14	2,10 mm	150 mm	20 mm
603 132	14	2,10 mm	200 mm	20 mm

Passende Schleuse
Fitting sheath

MRT Biopsie-Schleuse / MRI Biopsy Sheath			
Artikelnummer Order no.	Gauge	Durchmesser Diameter	Länge Length
601 408	16	1,60 mm	50 mm
601 410	16	1,60 mm	100 mm
601 412	16	1,60 mm	150 mm
601 418	15	1,90 mm	50 mm
601 420	15	1,90 mm	100 mm
601 422	15	1,90 mm	150 mm
601 428	13	2,40 mm	50 mm
601 430	13	2,40 mm	100 mm
601 432	13	2,40 mm	150 mm



SOMATEX[®]
Medical Technologies GmbH

Rheinstraße 7 d
D-14513 Teltow
Germany

Tel. +49 (0) 3328 30 76-0
Fax +49 (0) 3328 30 76-99
info@somatex.com
www.somatex.com